

Human IL-4 ELISA KIT

Catalog No. CDK029B

Quantity: 2 x 96 tests

PRODUCT SPECIFICATIONS :

Specificity : Recognizes both natural and recombinant human IL-4

Range : 1.1 pg / ml - 35 pg / ml

Sensitivity : < 0.5 pg/ml

Incubation : 3 h 50 min

Sample Types : Serum
Plasma
Cell culture supernatant

Cross Reaction : No cross reactivity with other human cytokines

Kit Content : Pre-coated 12 strip plate, biotinylated secondary antibody, standards, controls, buffers, Streptavidin-HRP, TMB Substrate, Stop Reagent.

1. INTENDED USE

The IL-4 ELISA kit is a solid phase sandwich ELISA for the *in-vitro* qualitative and quantitative determination of IL-4 in supernatants, buffered solutions or serum and plasma samples. This assay will recognize both natural and recombinant human IL-4.

This kit has been configured for research use only. Not suitable for use in therapeutic procedures.

2. PRINCIPLE OF THE METHOD

A capture Antibody highly specific for IL-4 has been coated to the wells of the microtiter strip plate provided during manufacture. Binding of IL-4 in samples and known standards to the capture antibodies is completed and then any excess unbound analyte is removed. During the next incubation period the binding of the biotinylated anti-IL-4 secondary antibody to the analyte occurs. Any excess unbound secondary antibody is then removed. The HRP conjugate solution is then added to every well including the zero wells, following incubation excess conjugate is removed by careful washing. A chromogen substrate is added to the wells resulting in the progressive development of a blue colored complex with the conjugate. The color development is then stopped by the addition of acid turning the resultant final product yellow. The intensity of the produced colored complex is directly proportional to the concentration of IL-4 present in the samples and standards. The absorbance of the color complex is then measured and the generated OD values for each standard are plotted against expected concentration forming a standard curve. This standard curve can then be used to accurately determine the concentration of IL-4 in any sample tested.



Cell Sciences[®]
480 Neponset Street
Bldg 12A
Canton, MA 02021

Toll Free: 888 769-1246
Phone: 781 828-0610
Fax: 781 828-0542

E-mail: info@cellsciences.com
Web Site: www.cellsciences.com

3. REAGENTS PROVIDED AND RECONSTITUTION

Reagents (Store @ 2-8°C)	Quantity 2 x 96 well kit	Reconstitution
P: 96 well microtiter strip plate	2	Ready to use (Pre-coated).
A: Plastic plate covers	4	n/a
B: Standard: 35 pg/ml	4	Reconstitute as directed on the vial (see Assay preparation, section 8).
C: Standard Diluent (Buffer)	1 (25 ml)	10x Concentrate, dilute in distilled water (see Assay preparation, section 8).
D: Standard Diluent (Human Serum)	2 (7.0 ml)	Ready to use.
E: Control	4	Reconstitute as directed on the vial (see Assay preparation, section 8).
F: Biotinylated anti-IL-4	2 (0.4 ml)	Dilute in Biotinylated Antibody Diluent (see Assay preparation, section 8).
G: Biotinylated Antibody Diluent	1 (13 ml)	Ready to use.
H: Streptavidin-HRP	4 (5 µl)	Add 0.5 ml of HRP Diluent prior to use (see Assay preparation, section 8).
I: HRP Diluent	1 (23 ml)	Ready to use.
J: Wash Buffer	2 (10 ml)	200x Concentrate dilute in distilled water (see Assay preparation, section 8).
K: TMB Substrate	1 (24 ml)	Ready to use.
L: H ₂ SO ₄ stop reagent	2 (11 ml)	Ready to use.

4. MATERIAL REQUIRED BUT NOT PROVIDED

- Microtiter plate reader fitted with appropriate filters (450 nm required with optional 630 nm reference filter)
- Microplate washer or wash bottle
- 10, 50, 100, 200 and 1,000 µl adjustable single channel micropipettes with disposable tips
- 50-300 µl multi-channel micropipette with disposable tips
- Multichannel micropipette reagent reservoirs
- Distilled water
- Vortex mixer
- Miscellaneous laboratory plastic and/or glass, if possible sterile

5. STORAGE INSTRUCTIONS

Store kit reagents between 2 and 8°C. Immediately after use, remaining reagents should be returned to cold storage (2-8°C). The expiration date of the kit and reagents is stated on box front labels. The expiration of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.



Cell Sciences®
480 Neponset Street
Bldg 12A
Canton, MA 02021

Toll Free: 888 769-1246
Phone: 781 828-0610
Fax: 781 828-0542

E-mail: info@cellsciences.com
Web Site: www.cellsciences.com

6. SPECIMEN COLLECTION, PROCESSING & STORAGE

Cell culture supernatants, human serum, plasma or other biological samples are suitable for use in the assay. Remove serum from the clot or red cells, respectively, as soon as possible after clotting and separation.

Cell culture supernatants: Remove particulates and aggregates by spinning at approximately 1000 x g for 10 min.

Serum: Use pyrogen/endotoxin free collecting tubes. Serum should be removed rapidly and carefully from the red cells after clotting. Following that, after clotting, centrifuge at approximately 1000 x g for 10 min and remove serum.

Plasma: EDTA, citrate, or heparin collected plasma can be assayed. Spin samples at 1000 x g for 30 min to remove particulates. Harvest plasma.

Storage: If not analyzed shortly after collection, samples should be aliquoted (250-500 µl) to avoid repeated freeze-thaw cycles and stored frozen at -80°C. Avoid multiple freeze-thaw cycles of frozen specimens.

Recommendation: Do not thaw by heating at 37°C or 56°C. Thaw at room temperature and make sure that the sample is completely thawed and homogeneous before use. When possible, avoid use of badly hemolyzed or lipemic sera. If large amounts of particles are present these should be removed prior to use by centrifugation or filtration.

7. SAFETY AND PRECAUTIONS FOR USE

- Handling of reagents, serum or plasma specimens should be in accordance with local safety procedures, e.g. CDC/NIH Health manual: "Biosafety in Microbiological and Biomedical Laboratories" 1984.
- The human serum included in this kit have been tested and found non-reactive for HBsAg, anti HIV1 & 2 and anti-VHC. Nevertheless, no known method can offer complete assurance that human blood derivatives will not transmit hepatitis, AIDS or other infections. Therefore handling of reagents, serum or plasma specimens should be in accordance with local safety procedures.
- Laboratory gloves should be worn at all times.
- Avoid any skin contact with H₂SO₄ and TMB substrate. In case of contact, wash thoroughly with water.
- Do not eat, drink, smoke or apply cosmetics where kit reagents are used.
- Do not pipette by mouth.
- When not in use, kit components should be stored refrigerated or frozen as indicated on vials or bottles labels.
- All reagents should be warmed to room temperature before use. Lyophilized standards should be discarded after use.
- Once the desired number of strips has been removed, immediately reseal the bag to protect the remaining strips from deterioration.
- Cover or cap all reagents when not in use.
- Do not mix or interchange reagents between different lots.
- Do not use reagents beyond the expiration date of the kit.
- Use a clean disposable plastic pipette tip for each reagent, standard, or specimen addition in order to avoid cross contamination. For the dispensing of H₂SO₄ and substrate solution, avoid pipettes with metal parts.
- Use a clean plastic container to prepare the washing solution.
- Thoroughly mix the reagents and samples before use by agitation or swirling.
- All residual washing liquid must be drained from the wells by efficient aspiration or by decantation followed by tapping the plate forcefully on absorbent paper. Never insert absorbent paper directly into the wells.



Cell Sciences[®]
480 Neponset Street
Bldg 12A
Canton, MA 02021

Toll Free: 888 769-1246
Phone: 781 828-0610
Fax: 781 828-0542

E-mail: info@cellsciences.com
Web Site: www.cellsciences.com

- The TMB solution is light sensitive. Avoid prolonged exposure to light. Also, avoid contact of the TMB solution with metal to prevent color development. **Warning TMB is toxic avoid direct contact with hands. Dispose of properly.**
- If a dark blue color develops within a few minutes after preparation, this indicates that the TMB solution has been contaminated and must be discarded. Read absorbances immediately after completion of the assay.
- When pipetting reagents, maintain a consistent order of addition from well-to-well. This will ensure equal incubation times for all wells.
- Follow incubation times described in the assay procedure.
- Dispense the TMB solution immediately after washing microtiter plate.

8. ASSAY PREPARATION

Bring all reagents to room temperature before use.

8.1 Assay Design

Determine the number of microwell strips required to test the desired number of samples plus the appropriate number of wells needed for running zeros and standards. Each sample, standard, and zero should be tested **in duplicate**. Remove sufficient Microwell Strips for testing from the pouch immediately prior to use. Return any wells not required for this assay with desiccant to the pouch. Seal tightly and return to 2-8 °C storage.

Example plate layout (example shown for a 6 point standard curve)

	Standards		Sample Wells									
	1	2	3	4	5	6	7	8	9	10	11	12
A	35	35										
B	17.5	17.5										
C	8.75	8.75										
D	4.4	4.4										
E	2.1	2.1										
F	1	1										
G	zero	zero										
H	control	control										

All remaining empty wells can be used to test samples in duplicate.

8.2 Preparation of Wash Buffer

Dilute the (200x) wash buffer concentrate 200 fold with distilled water to give a 1x working solution. Pour entire contents (10 ml) of the Washing Buffer Concentrate into a clean 2,000 ml graduated cylinder. Bring final volume to 2,000 ml with glass-distilled or deionized water. Mix gently to avoid foaming. Transfer to a clean wash bottle and store at 2°-8°C.

8.3 Preparation of Standard Diluent Buffer

Add the contents of the vial (10x concentrate) to 225 ml of distilled water before use. This solution can be stored at 2-8 °C for up to 1 week.



Cell Sciences[®]
480 Neponset Street
Bldg 12A
Canton, MA 02021

Toll Free: 888 769-1246
Phone: 781 828-0610
Fax: 781 828-0542

E-mail: info@cellsciences.com
Web Site: www.cellsciences.com

8.4 Preparation of Standard

Depending on the type of samples you are assaying, the kit includes two standard diluents because biological fluids might contain proteases or cytokine-binding proteins that could modify the recognition of the cytokine you want to measure. You should reconstitute standard vials with the most appropriate standard diluent.

For **serum and plasma** samples: use Standard Diluent for human serum.

For **cell culture supernatants**: use Standard Diluent Buffer.

Standard vials must be reconstituted with the volume of standard diluent shown on the vial immediately prior to use. This reconstitution gives a stock solution of 35 pg/ml of IL-4. Mix the reconstituted standard gently by repeated aspirations/ejections. Serial dilutions of the standard are made directly in the assay plate to provide the concentration range from 35 to 1 pg/ml. A fresh standard curve should be produced for each new assay.

- Immediately after reconstitution, add 200 μ l of the reconstituted standard to wells A1 and A2, which provides the highest concentration standard at 35 pg/ml.
- Add 100 μ l of appropriate Standard Diluent to the remaining standard wells B1 and B2 to F1 and F2.
- Transfer 100 μ l from wells A1 and A2 to B1 and B2. Mix the well contents by repeated aspirations and ejection taking care not to scratch the inner surface of the wells.
- Continue this 1:1 dilution using 100 μ l from wells B1 and B2 through to wells F1 and F2 providing a serial diluted standard curve ranging from 35 pg/ml to 1 pg/ml.
- Discard 100 μ l from the final wells of the standard curve (F1 and F2).

Alternatively these dilutions can be performed in separate clean tubes and immediately transferred directly into the relevant wells.

8.5. Preparation of Control

Freeze-dried control vials should also be reconstituted with the most appropriate standard diluent to your samples.

For **serum and plasma** samples: use Standard diluent - human serum

For **cells culture supernatants**: use Standard diluent buffer.

The Control must be reconstituted with the volume of standard buffer diluent indicated on the vial. Reconstitution of the freeze-dried material with the indicated volume will give a solution for which the IL-4 concentration is stated on the vial. Do not store after use.



8.6. Preparation of Biotinylated anti-IL-4

It is recommended this reagent is prepared immediately before use. Dilute the biotinylated anti-IL-4 with the biotinylated antibody diluent in an appropriate clean glass vial using volumes appropriate to the number of required wells. Please see example volumes below:

Number of wells required	Biotinylated Antibody (μ l)	Biotinylated Antibody Diluent (μ l)
16	40	1060
24	60	1590
32	80	2120
48	120	3180
96	240	6360

8.7. Preparation of Streptavidin-HRP

It is recommended to centrifuge vial for a few seconds in a microcentrifuge to collect all the volume at the bottom.

Dilute the 5 μ l vial with 0.5 ml of HRP diluent **immediately before use**. Do not keep this diluted vial for future experiments. Further dilute the HRP solution to volumes appropriate for the number of required wells in a clean glass vial. Please see example volumes below:

Number of wells required	Streptavidin-HRP (μ l)	Streptavidin-HRP Diluent (μ l)
16	30	2
24	45	3
32	60	4
48	75	5
96	150	10



9. METHOD

We strongly recommend that every vial is mixed thoroughly without foaming prior to use.

Prepare all reagents as shown in section 8.

Note: Final preparation of Biotinylated anti-IL-4 (section 8.6) and Streptavidin-HRP (section 8.7) should occur immediately before use.

Assay Step		Details
1.	Addition	Prepare Standard curve as shown in section 8.4 above
2.	Addition	Add 100 µl of each standard, sample, control, and zero (appropriate standard diluent) in duplicate to the appropriate number of wells.
3.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for 2 hours .
4.	Wash	Remove the cover and wash the plate as follows: a) Aspirate the liquid from each well. b) Dispense 0.3 ml of 1x washing solution into each well. c) Aspirate the contents of each well. d) Repeat step b and c another two times.
5.	Addition	Add 50 µl of diluted biotinylated anti-IL-4 to all wells.
6.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for 1 hour .
7.	Wash	Repeat wash step 4.
8.	Addition	Add 100 µl of Streptavidin-HRP solution into all wells.
9.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25°C) for 30 min
10.	Wash	Repeat wash step 4.
11.	Addition	Add 100µl of ready-to-use TMB Substrate Solution into all wells.
12.	Incubation	Incubate in the dark for 12-15 minutes* at room temperature. Avoid direct exposure to light by wrapping the plate in aluminum foil.
13.	Addition	Add 100µl of H2SO4:Stop Reagent into all wells.
<p>Read the absorbance value of each well (immediately after step 11) on a spectrophotometer using 450 nm as the primary wave length and optionally 620 nm as the reference wave length (610 nm to 650 nm is acceptable).</p>		

**Incubation time of the substrate solution is usually determined by the ELISA reader performance. Many ELISA readers only record absorbance up to 2.0 O.D. Therefore the color development within individual microwells must be observed by the analyst, and the substrate reaction stopped before positive wells are no longer within recordable range.*



Cell Sciences®
 480 Neponset Street
 Bldg 12A
 Canton, MA 02021

Toll Free: 888 769-1246
 Phone: 781 828-0610
 Fax: 781 828-0542

E-mail: info@cellsciences.com
 Web Site: www.cellsciences.com

10. DATA ANALYSIS


Calculate the average absorbance values for each set of duplicate standards and samples. Ideally duplicates should be within 20% of the mean.

Generate a linear standard curve by plotting the average absorbance of each standard on the vertical axis versus the corresponding Human IL-4 standard concentration on the horizontal axis.

The amount of IL-4 in each sample is determined by extrapolating OD values against IL-4 standard concentrations using the standard curve.

Example IL-4 Standard Curve

Standard	IL-4 Concentration (pg/ml)	OD (450 nm) mean	CV (%)
1	35	2.559	2.4
2	17.5	1.347	2.7
3	8.75	0.672	4.0
4	4.4	0.346	0.2
5	2.1	0.204	2.6
6	1	0.139	4.6
Zero	0	0.066	13.5

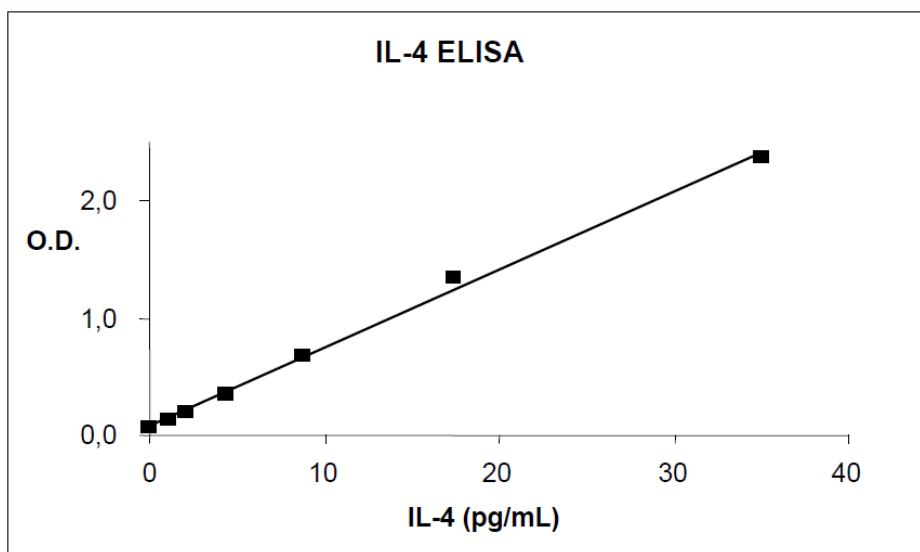

Note: Curve shown above should not be used to determine results. Every laboratory must produce a standard curve for each set of microwell strips assayed.



Cell Sciences[®]
480 Neponset Street
Bldg 12A
Canton, MA 02021

Toll Free: 888 769-1246
Phone: 781 828-0610
Fax: 781 828-0542

E-mail: info@cellsciences.com
Web Site: www.cellsciences.com

11. ASSAY LIMITATIONS

Do not extrapolate the standard curve beyond the maximum standard curve point. The dose-response is non-linear in this region and good accuracy is difficult to obtain. Concentrated samples above the maximum standard concentration must be diluted with Standard diluent or with your own sample buffer to produce an OD value within the range of the standard curve. Following analysis of such samples, always multiply results by the appropriate dilution factor to produce the actual final concentration.

The influence of various drugs on end results has not been investigated. Bacterial or fungal contamination and laboratory cross-contamination may also cause irregular results.

Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Completely empty wells before dispensing fresh Washing Buffer. Fill wells with Washing Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.

Disposable pipette tips, flasks or glassware are preferred. Reusable glassware must be washed and thoroughly rinsed of all detergents before use.

As with most biological assays, conditions may vary from assay to assay. Therefore, **a fresh standard curve must be prepared and run for every assay.**

12. PERFORMANCE CHARACTERISTICS

12.1 Sensitivity

The minimum detectable dose of IL-4 was determined to be **0.7 pg/ml**.

This was determined by adding 3 standard deviations to the mean optical density obtained when the zero standard was assayed 35 times.

12.2 Specificity

The assay recognizes both natural and recombinant human IL-4. To define the specificity of this ELISA several proteins were tested for cross reactivity. There was no cross reactivity observed for any protein tested (IL-1 α , IL-1 β , IL-10 IL-12, IFN γ , IL-2, IL-6, TNF α , IL-8 and IL-13).

12.3 Precision

Intra-assay

Reproducibility within the assay was evaluated in three independent experiments. Each assay was carried out with 6 replicates (3 duplicates) of 3 spiked serum or human pooled serum samples and 2 supernatants containing different concentrations of IL-4. **The overall intra-assay coefficient of variation has been calculated to be 3.5%.**



Session	Sample	Assay 1 IL-4 pg/ml	Assay 2 IL-4 pg/ml	Assay 3 IL-4 pg/ml	Mean IL-4 pg/ml	SD	CV %
1	1	30	28	29	29	1.0	3.4
	2	21	20	20	20	0.6	2.8
	3	10	11	11	11	0.6	5.4
	4	30	29	29	29	0.6	1.9
	5	22	21	22	22	0.6	2.6
2	1	23	23	23	23	0.4	1.6
	2	15	16	16	16	0.2	0.9
	3	10	12	10	11	1.3	11.8
	4	28	28	27	28	0.4	1.3
	5	22	22	22	22	0.4	1.8
3	1	26	26	24	25	1.2	4.6
	2	18	18	18	18	0	0
	3	11	11	13	12	1.2	10
	4	32	32	31	32	0.6	1.8
	5	24	24	25	24	0.6	2.4

Inter-assay

Assay to assay reproducibility within one laboratory will be evaluated in three independent experiments by two technicians. Each assay will be carry out with 6 replicates of 3 spiked serum human pooled serum samples and 2 supernatants containing different concentration of IL-4. **The calculated overall coefficient of variation was 7.3%.**

Technician	Session	Sample 1 IL-4 pg/ml	Sample 2 IL-4 pg/ml	Sample 3 IL-4 pg/ml	Sample 4 IL-4 pg/ml	Sample 5 IL-4 pg/ml
A	1	30	21	10	30	22
		28	20	11	29	21
		29	20	11	29	22
	2	23	16	9	28	22
		23	16	12	28	22
		23	16	10	27	22
	3	26	18	11	32	24
		26	18	11	32	24
		24	18	13	31	25
B	1	24	17	13	28	21
		23	15	11	28	22
		22	15	10	27	11
	2	27	18	14	31	23
		25	18	13	30	24
		24	18	12	29	24
	3	26	17	12	27	21
		24	17	11	27	21
		23	16	11	27	23
Mean IL-4 pg/ml		24	17	11	29	23
SD		1	1	1	2	1
CV %		6.2	6.6	11.1	6.6	6.3



12.4 Dilution Parallelism

Four human pooled serum samples with different levels of IL-4 were analyzed at different serial two fold dilutions (1:2 to 1:8) with two replicates each. Recoveries ranged from 75 to 121% with an overall **mean recovery of 100%**.

12.5 Spike Recovery

The spike recovery was evaluated by spiking 3 concentrations of IL-4 in human serum in 3 separate experiments. Recoveries ranged from 100 to 113% with an overall **mean recovery of 107%**.

12.6. Stability

Storage Stability

Aliquots of spiked serum samples were stored at -20°C , $2-8^{\circ}\text{C}$, room temperature (RT) and at 37°C and the IL-4 level determined after 24 h. There was no significant loss of IL-4 reactivity during storage at RT, $2-8^{\circ}\text{C}$ or 37°C .

Freeze-thaw Stability

Aliquots of spiked serum were stored frozen at -20°C and thawed, refrozen, and re-thawed up to 5 times and the IL-4 level was determined. There was no significant loss of IL-4 after 5 cycles of freezing and thawing.

12.7. Expected serum values

A panel of 16 human sera was tested for IL-4. All were below the detection level of 0.7 pg/ml.


12.8. Standard Calibration

This immunoassay is calibrated against the International Reference Standard NIBSC 88/656. NIBSC 88/656 is quantitated in International Units (IU) and equivalence in ng/ml is indicated. It has been calculated that 1 ng correspond to 1 ng calibration IL-4.



13. ASSAY SUMMARY

Total procedure length: 3 h 45 min

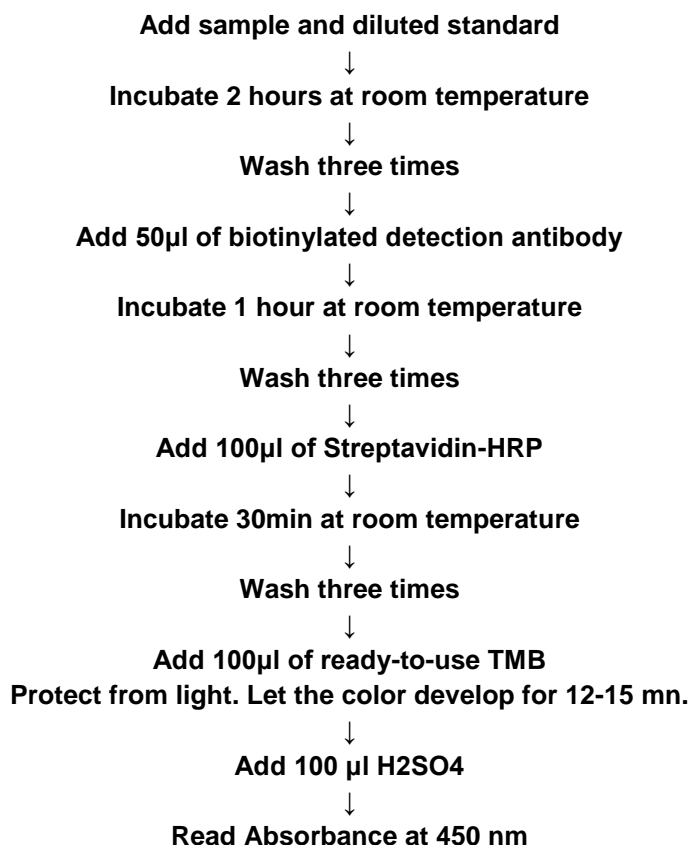

NOT FOR HUMAN USE. FOR RESEARCH ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE



Cell Sciences[®]
480 Neponset Street
Bldg 12A
Canton, MA 02021

Toll Free: 888 769-1246
Phone: 781 828-0610
Fax: 781 828-0542

E-mail: info@cellsciences.com
Web Site: www.cellsciences.com