

## Human CD102 / ICAM-2 ELISA Kit

**Catalog No.** CDK020B

**Quantity:** 2 x 96 tests

### PRODUCT SPECIFICATIONS :

**Specificity :** Recognizes both natural and recombinant human soluble ICAM-2

**Range :** 0.75 U / ml – 24 U / ml

**Sensitivity :** 0.2 U / ml

**Incubation :** 1 h 40 min

**Sample Types :** Serum  
Plasma  
Cell culture supernatant

**Cross Reaction :** No cross reactivity with other human adhesion molecules

**Kit Content :** Pre-coated 12 strip plates, biotinylated secondary antibody, standards, controls (when available), buffers, Streptavidin-HRP, TMB, Stop Reagent.

### 1. INTENDED USE

The Cell Sciences Inc.® Human sICAM-2/CD102 ELISA is to be used for the *in-vitro* quantitative determination of soluble Inter Cellular Adhesion Molecule -2 (sICAM-2) in human serum, plasma, buffered solutions or cell culture medium. The assay will recognize both natural and recombinant human sICAM-2/CD102.

**This kit has been configured for research use only.**

### 2. PRINCIPLE OF THE METHOD

The sICAM-2/CD102 Kit is a solid phase sandwich Enzyme Linked-Immuno- Sorbent Assay (ELISA). A monoclonal antibody specific for ICAM-2 has been coated onto the wells of the microtiter strips provided. Samples, including standards of known sICAM-2 concentrations, controls and unknowns are pipetted into these wells.

During the first incubation, the sICAM-2 antigen and a biotinylated monoclonal antibody specific for ICAM-2 are simultaneously incubated.

After washing, the enzyme (streptavidin-peroxydase) is added. After incubation and washing to remove the unbound enzyme, a substrate solution is added to induce a colored reaction product. The intensity of this colored product is directly proportional to the concentration of sICAM-2 present in the samples.



**Cell Sciences, Inc.®**  
65 Parker Street  
Unit 11  
Newburyport, MA 01950

Toll Free: 888-769-1246  
Phone: 978-572-1070  
Fax: 978-992-0298

E-mail: [info@cellsciences.com](mailto:info@cellsciences.com)  
Web Site: [www.cellsciences.com](http://www.cellsciences.com)

### 3. REAGENTS PROVIDED AND RECONSTITUTION

Reagents (Store @ 2-8°C)	Quantity 2 x 96 well kit	Reconstitution
P: 96 well microtiter strip plate	2	Ready to use. (Pre-coated)
A: Plastic plate covers	4	n/a
B: Standard: 24 U/ml	4 vials	Reconstitute as directed on the vial. (see Assay preparation, section 8)
C: Control	4 vials	Reconstitute as directed on the vial. (see Assay preparation, section 8)
D: Standard Diluent (Buffer)	1 vial (25 ml)	10x Concentrate, dilute in distilled water (see reagent preparation, section 8).
E: Biotinylated anti-ICAM-2	2 vial ( 0.4 ml)	Dilute in biotinylated antibody diluent (see reagent preparation, section 8).
F: Biotinylated Antibody diluent	1 vial (13 ml)	Ready to use.
G: Streptavidin-HRP	4 vials (5 µl)	Add 0.5 ml of HRP diluent prior to use. (see Assay preparation, section 8)
H: HRP Diluent	1 vial (23 ml)	Ready to use.
I: Wash Buffer	2 vials (10 ml)	200x Concentrate, dilute in distilled water. (see Assay preparation, section 8)
J: TMB Substrate	1 vial (24 ml)	Ready to use.
K: H <sub>2</sub> SO <sub>4</sub> stop reagent	2 vials (11 ml)	Ready to use.

### 4. MATERIAL REQUIRED BUT NOT PROVIDED

- Microtiter plate reader fitted with appropriate filters (450 nm required with optional 620 nm reference filter)
- Microplate washer or wash bottle
- 10, 50, 100, 200 and 1,000 µl adjustable single channel micropipettes with disposable tips
- 50-300 µl multi-channel micropipette with disposable tips
- Multichannel micropipette reagent reservoirs
- Distilled water
- Vortex mixer
- Miscellaneous laboratory plastic and/or glass, if possible sterile



## 5. STORAGE INSTRUCTIONS

Store the kit reagents between 2 and 8 °C. Immediately after use, remaining reagents should be returned to cold storage (2-8 °C). The expiration date of the kit and reagents is stated on box front labels. The expiration of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

**Wash Buffer 1X:** Once prepared, store at 2-8 °C for up to 1 week.

**Standard Diluent Buffer 1X:** Once prepared, store at 2-8 °C for up to 1 week.

**Reconstituted Standards / Controls:** Once prepared, use immediately and do not store.

**Diluted Biotinylated Secondary Antibody:** Once prepared, use immediately and do not store.

**Diluted Streptavidin-HRP:** Once prepared, use immediately and do not store.

## 6. SPECIMEN COLLECTION, PROCESSING & STORAGE

Cell culture supernatants, serum, plasma or other biological samples are suitable for use in the assay. Remove serum from the clot or red cells respectively, as soon as possible after clotting and separation.

**Cell culture supernatants:** Remove particulates and aggregates by spinning at approximately 1000 x g for 10 min.

**Serum:** Avoid any unintentional stimulation of the cells by the procedure. Use pyrogen/endotoxin free collecting tubes. Serum should be removed rapidly and carefully from the red cells after clotting. After clotting, centrifuge at approximately 1000 x g for 10 min and remove serum.

**Before testing, human serum samples must be diluted 20 times in standard buffer diluent.**

**Plasma:** EDTA, citrate and heparin plasma can be assayed. Spin samples at 1000 x g for 30 min. to remove particulates. Harvest plasma.

**Storage:** If not analyzed shortly after collection, samples should be aliquoted (250-500 µl) to avoid repeated freeze-thaw cycles and stored frozen at -80 °C. Avoid multiple freeze-thaw cycles of frozen specimens.

**Recommendation:** Do not thaw by heating at 37 °C or 56 °C. Thaw at room temperature, and make sure that the sample is completely thawed and homogeneous before use. When possible, avoid use of badly hemolyzed or lipemic sera. If large amounts of particles are present, these should be removed prior to use by centrifugation or filtration.



## 7. SAFETY AND PRECAUTIONS FOR USE

- Handling of reagents, serum, or plasma specimens should be in accordance with local safety procedures, e.g. CDC/NIH Health manual: "Biosafety in Microbiological and Biomedical Laboratories" 1984.
- The human serum included in this kit has been tested, and found non-reactive for HBsAg, anti HIV1 & 2 and anti VHC. Nevertheless, no known method can offer complete assurance that human blood derivatives will not transmit hepatitis, AIDS or other infections. Therefore, handling of reagents, serum or plasma specimens should be in accordance with local safety procedures.
- Laboratory gloves should be worn at all times.
- Avoid any skin contact with H<sub>2</sub>SO<sub>4</sub> and TMB. In case of contact, wash thoroughly with water.
- Do not eat, drink, smoke or apply cosmetics where kit reagents are used.
- Do not pipette by mouth.
- When not in use, kit components should be stored refrigerated or frozen, as indicated on vial or bottle labels.
- All reagents should be warmed to room temperature before use. Lyophilized standards should be discarded after use.
- Once the desired number of strips has been removed, immediately reseal the bag to protect the remaining strips from deterioration.
- Cover or cap all reagents when not in use.
- Do not mix or interchange reagents between different lots.
- Do not use reagents beyond the expiration date of the kit.
- Use a clean, disposable, plastic pipette tip for each reagent, standard, or specimen addition in order to avoid cross contamination for the dispensing of H<sub>2</sub>SO<sub>4</sub> and substrate solution. Avoid pipettes with metal parts.
- Use a clean plastic container to prepare the washing solution.
- Thoroughly mix the reagents and samples before use by agitation or swirling.
- All residual washing liquid must be drained from the wells by efficient aspiration, or by decantation, followed by tapping the plate forcefully on absorbent paper. Never insert absorbent paper directly into the wells.
- The TMB solution is light sensitive. Avoid prolonged exposure to light. Also, avoid contact of the TMB solution with metal to prevent color development. **Warning: TMB is toxic. Avoid direct contact with hands. Dispose of properly.**
- If a dark blue color develops within a few minutes after preparation, this indicates that the TMB solution has been contaminated and must be discarded. Read absorbances within 1 hour after completion of the assay.
- When pipetting reagents, maintain a consistent order of addition from well-to-well. This will ensure equal incubation times for all wells.
- Follow incubation times described in the assay procedure.
- Dispense the TMB solution immediately after washing of the microtiter plate.



## 8. ASSAY PREPARATION

Bring all reagents to room temperature before use.

### 8.1 Assay Design

Determine the number of microwell strips required to test the desired number of samples, plus appropriate number of wells needed for running zeros and standards. Each sample, standard and zero should be tested **in duplicate**. Remove sufficient Microwell Strips for testing from the pouch immediately prior to use. Return any wells not required for this assay with desiccant to the pouch. Seal tightly and return to 2-8°C storage.

**Example plate layout** (example shown for a 6 point standard curve)

	Standards (U/mL)		Sample Wells									
	1	2	3	4	5	6	7	8	9	10	11	12
A	24	24										
B	12	12										
C	6	6										
D	3	3										
E	1.5	1.5										
F	0.75	0.75										
G	zero	zero										
H	ctrl	ctrl										

*All remaining empty wells can be used to test samples in duplicate*

### 8.2 Preparation of Wash Buffer

Dilute the (200x) wash buffer concentrate 200 fold with distilled water to give a 1x working solution. Pour entire contents (10 ml) of the Washing Buffer Concentrate into a clean 2,000 ml graduated cylinder. Bring final volume to 2,000 ml with glass-distilled or deionized water. Mix gently to avoid foaming. Transfer to a clean wash bottle and store at 2°-8°C for up to 1 week.

### 8.3 Preparation of Standard Diluent Buffer

Add the contents of the vial (10x concentrate) to 225 ml of distilled water before use.

This solution can be stored at 2-8°C for up to 1 week.

### 8.4 Preparation of Standard

Standard vials must be reconstituted with the volume of standard diluent shown on the vial immediately prior to use. This reconstitution gives a stock solution of 10 U/ml of sICAM-2. **Mix the reconstituted standard gently by inversion only.** Serial dilutions of the standard are made directly in the assay plate to provide the concentration range from 24 to 0.75 U/ml. A fresh standard curve should be produced for each new assay.

- Immediately after reconstitution add 200µl of the reconstituted standard to wells A1 and A2, which provides the highest concentration standard at 24 U/ml



- Add 100µl of Standard Diluent to the remaining standard wells B1 and B2 to F1 and F2
- Transfer 100µl from wells A1 and A2 to B1 and B2. Mix the well contents by repeated aspirations and ejections taking care not to scratch the inner surface of the wells
- Continue this 1:1 dilution using 100µl from wells B1 and B2 through to wells F1 and F2 providing a serial diluted standard curve ranging from 24 U/ml to 0.75 U/ml
- Discard 100µl from the final wells of the standard curve (F1 and F2)

Alternatively these dilutions can be performed in separate clean tubes and immediately transferred directly into the relevant wells.

## 8.5 Preparation of Controls

Controls must be reconstituted with the volume of standard buffer diluent indicated on the vial. Reconstitution of the freeze-dried material with the indicated volume will give a solution for which the sICAM-2 concentration is stated on the vial. Allow control to stand for 5 minutes with gentle swirling prior to distributing in control wells. Do not store after use.

## 8.6 Preparation of Samples

Before testing, human serum or samples must be diluted 20 times in standard diluent buffer.

## 8.7 Preparation of Biotinylated anti-sICAM-2

It is recommended this reagent is prepared immediately before use. Dilute the biotinylated anti-sICAM-2 with the biotinylated antibody diluent in an appropriate clean glass vial using volumes appropriate to the number of required wells. Please see example volumes below:

Number of wells required	Biotinylated Antibody (µl)	Biotinylated Antibody Diluent (µl)
16	40	1060
24	60	1590
32	80	2120
48	120	3180
96	240	6360

## 8.8 Preparation of Streptavidin-HRP

It is recommended to centrifuge vial for a few seconds in a microcentrifuge to collect all the volume at the bottom.

Dilute the 5 µl vial with 0.5 ml of HRP diluent **immediately before use**. Do-not keep this diluted vial for future experiments. Further dilute the HRP solution to volumes appropriate for the number of required wells in a clean glass vial. Please see example volumes below:



Number of wells required	Streptavidin-HRP ( $\mu$ l)	Streptavidin-HRP Diluent (ml)
16	30	2
24	45	3
32	60	4
48	75	5
96	150	10

## 9. METHOD

We strongly recommend that every vial is mixed thoroughly without foaming prior to use except the standard vial which must be mixed gently by inversion only.

Prepare all reagents as shown in section 8. Note: Final preparation of Biotinylated anti-sICAM-2 (section 8.7) and Streptavidin-HRP (section 8.8) should occur immediately before use.

Assay Step		Details
1.	Preparation	<b>Prepare Standard Curve</b> as shown in 8.4
2.	Addition	Add 100 $\mu$ l of each <b>standard, sample, control and zero</b> (standard diluent) in duplicate to appropriate number of wells
3.	Addition	Add 50 $\mu$ l of diluted <b>Biotinylated anti-sICAM-2</b> to all wells
4.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25 °C) for <b>1 hour</b> .
5.	Wash	Remove the cover and wash the plate as follows: a) Aspirate the liquid from each well. b) Dispense 0.3 ml of <b>1x washing solution</b> into each well. c) Aspirate the contents of each well. d) Repeat step b and c another two times.
6.	Addition	Add 100 $\mu$ l of <b>Streptavidin-HRP</b> into all wells.
7.	Incubation	Cover with a plastic plate cover and incubate at room temperature (18 to 25 °C) for <b>30 minutes</b> .
8.	Wash	Repeat wash as in Step 5.
9.	Addition	Add 100 $\mu$ l of ready-to-use <b>TMB Substrate Solution</b> into all wells.
10.	Incubation	Incubate in the dark for <b>10-20 minutes*</b> at room temperature. Avoid direct exposure to light by wrapping the plate in aluminum foil.
11.	Addition	Add 100 $\mu$ l of <b>H<sub>2</sub>SO<sub>4</sub>: Stop Reagent</b> into all wells.



**Read the absorbance** value of each well (immediately after step 11) on a spectrophotometer, using 450 nm as the primary wavelength, and optionally 630 nm as the reference wave length (610 nm to 650 nm is acceptable).

*\*Incubation time of the substrate solution is usually determined by the ELISA reader performance. Many ELISA readers only record absorbance up to 2.0 O.D. Therefore the color development within individual microwells must be observed by the analyst, and the substrate reaction stopped before positive wells are no longer within recordable range.*

## 10. DATA ANALYSIS


Calculate the average absorbance values for each set of duplicate standards and samples. Ideally, duplicates should be within 20% of the mean.

Generate a linear standard curve by plotting the average absorbance of each standard on the vertical axis versus the corresponding Human sICAM-2 standard concentration on the horizontal axis.

The amount of sICAM-2 in each sample is determined by extrapolating O.D. values against sICAM-2 standard concentrations using the standard curve.

### Sample ICAM-2 Curve:

Standard	ICAM-2 Conc U/ml	OD (450nm) Mean	CV (%)
1	24	1.508	1.3
2	12	0.863	2.5
3	6	0.527	6.8
4	3	0.289	0.2
5	1.5	0.141	9.0
6	0.75	0.111	3.9
Zero	0	0.038	-

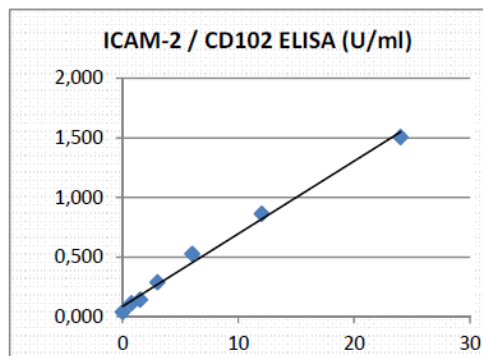

**Note:** The curve shown above should not be used to determine results. Every lab must produce a standard curve for each set of microwell strips assayed.

For serum and plasma samples which have been diluted according to the protocol (1:20), the calculated concentration should be multiplied by the dilution factor (x20).

## 11. ASSAY LIMITATIONS

Do not extrapolate the standard curve beyond the maximum standard curve point. The dose-response is non-linear in this region and good accuracy is difficult to obtain. Concentrated samples above the maximum standard concentration must be diluted with Standard diluent or with your own sample buffer to produce an OD value within the range of the standard curve. Following analysis of such samples always multiply results by the appropriate dilution factor to produce actual final concentration.

The influence of various drugs on end results has not been investigated. Bacterial or fungal contamination and laboratory cross-contamination may also cause irregular results.

Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Completely empty wells before dispensing fresh Washing Buffer, fill with Washing Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.

Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.

As with most biological assays conditions may vary from assay to assay therefore **a fresh standard curve must be prepared and run for every assay.**

## 12. PERFORMANCE CHARACTERISTICS

### 12.1 Sensitivity

The minimum detectable dose of sICAM-2 is less than **0.2 U/ml**.

This has been determined by adding 3 standard deviations to the mean optical density obtained when the zero standard was assayed 36 times.

### 12.2. Precision

Intra-Assay					Inter-Assay				
Sample	n	Mean (U/mL)	SD	CV%	Sample	n	Mean (U/mL)	SD	CV%
A	8	23.62	0.277	1.17	A	33	23.26	0.668	2.87
B	8	2.96	0.086	2.90	B	33	3.12	0.218	6.98

### 12.3. Normal Serum Values

The average concentration of sICAM-2 detected in 80 normal human serum is 220+/-74 U/ml ranging from 44 to 405 U/ml.



## 13. ASSAY SUMMARY

### Total procedure length : 1h 45mn

Add 100 µl of sample or diluted standard or control



Add 50µl of diluted biotinylated  
detection antibody to all wells



Incubate 1 hour at room temperature



Wash three times



Add 100 µl of streptavidin-HRP to all wells



Incubate 30 min at room temperature



Wash three times



Add 100 µl of ready-to-use TMB  
Protect from light. Let the color develop for 10-20 min.



Add 100 H<sub>2</sub>SO<sub>4</sub>



Read Absorbance at 450 nm

**NOT FOR HUMAN USE. FOR RESEARCH ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.**



**Cell Sciences, Inc.**<sup>®</sup>  
65 Parker Street  
Unit 11  
Newburyport, MA 01950

Toll Free: 888-769-1246  
Phone: 978-572-1070  
Fax: 978-992-0298

E-mail: [info@cellsciences.com](mailto:info@cellsciences.com)  
Web Site: [www.cellsciences.com](http://www.cellsciences.com)